Biochemistry lab

Experiment (8)

EXPIREMNT 8

ENZYMATIC COLORIMETRIC METHOD (END POINT)


Purpose: to determine the concentration of components inside the serum by determining its absorbance, using spectrophotometer.

- Serum: plasma without RBC.
- Serum components: 1-Triglycerides (lipids)
2-glucose (carbohydrates)
3-cholesterol (lipids)
4-albumin (proteins)

Purpose of the experiment: measuring those components' concentrations.

*we will measure the concentration by beer's law:

Remember: $A = (a.b).c$ $\longleftrightarrow A1/C1 = A2/C2$

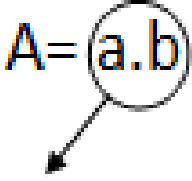Constant

A: Absorbance.

C: Concentration.

a: absorption coefficient.

b: cuvette's width.

- First, here are some abbreviations you should know:

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

1st experiment: measuring the concentration of triglycerides.



Triglycerides
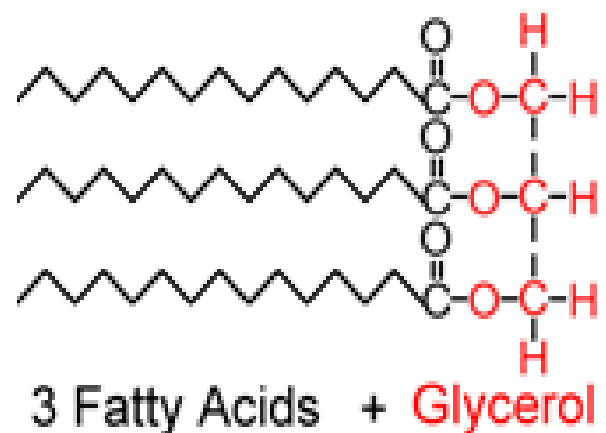
$$\lambda_{\max} = 500 \text{ nm}$$

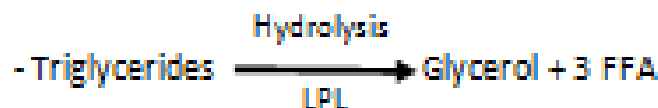Those should be remembered.

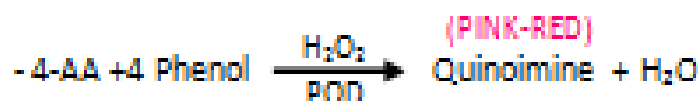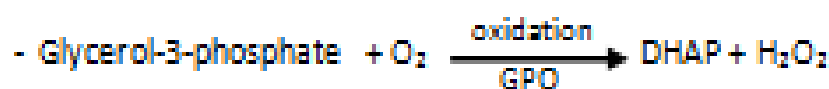
- as color's intensity increases, [H₂O₂] increases.
- [H₂O₂] describes the concentration of Triglyceride.
- Triglyceride's reagent consists of: LPL, GK, GPO, POD. When it reacts with triglyceride it gives us Quinoimine (pink – red)

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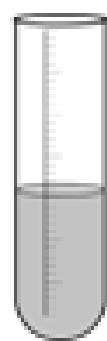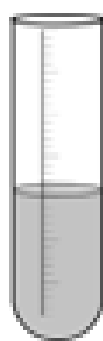
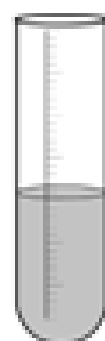
10 μ l

standard (200mg/dl)

Remember: blank solution is the solvent without any solute; we reset the device on it to neglect the solvent's absorption.

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Procedure:

- Prepare three test tubes then add 1ml of reagent to each one of them.
- Add 10 μ l of the sample (serum) to the second test tube.
- Add 10 μ l of the standard to the third test tube.
- Measure the absorbance by spectrophotometer at 37 degrees at 500nm wavelength.
- Find the concentration of the sample by beer's law.

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

$$C_{\text{standard}} = 200 \text{ mg/dl}$$

$$C_{\text{sample}} = (A_{\text{sample}} * C_{\text{standard}}) / A_{\text{standard}}$$

2nd experiment: Glucose

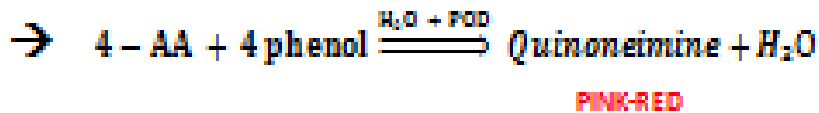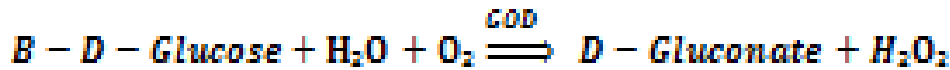Oxidation $\rightarrow H_2O_2$ appears $\rightarrow$ Quinoneimine appearance.

Glucose's reagent consists of: GOD, POD. It gives Quinoneimine (pink).

● Practical part:

- Add 3 test tubes each with 1 ml of reagent
- Add of serum to the test tube number 2
- Add 10 μ l of standard glucose to test tube number 3
- Because the first one is blank, contains the reagent only. Don't add anything else.
- Measure the absorbance at 37° C, $\lambda = 500\text{nm}$
- Then measure the concentration by beer's law.

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
$$C_{\text{standard}} = 100\text{mg/dl}$$

3rd Experiment: - Cholesterol

- Cholesterol is found in two shapes: Free (minor)

Esterfied (major)

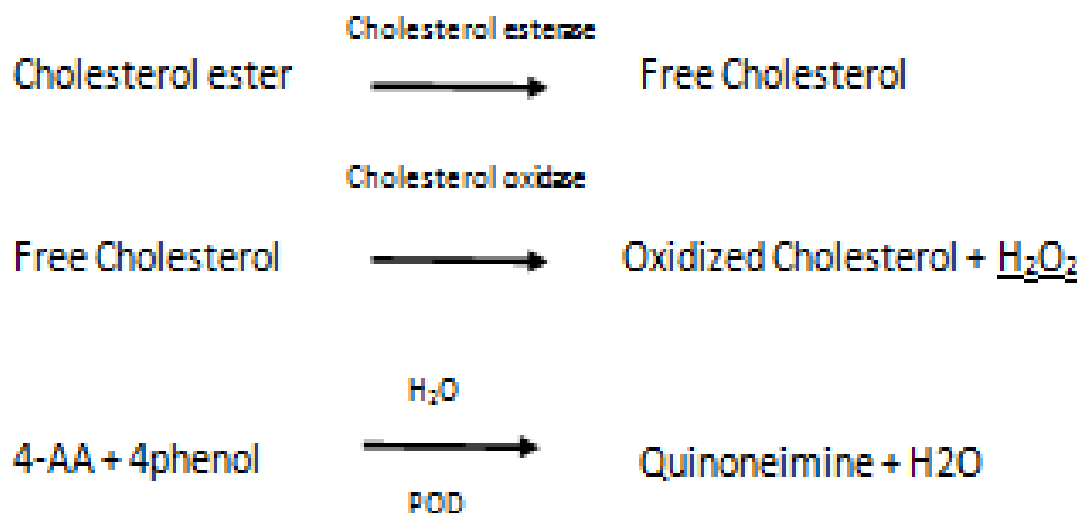- Cholesterol's agent consists of: Cholesterol esterase, POD

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Practical part:-

- Add 3 test tubes with 1mL of reagent on each one.
- Test tube 1 will be Blank, add 10 μ L of the sample (serum) to test tube 2, and 10 μ L of standard to test tube 3.
- Then measure the absorbance at 37° and at
- Then measure the concentration by Beer's Law.

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
$$C_{\text{standard}} = 200 \text{ mg/dl}$$

4th Experiment: - proteins

- Albumin is considered 60 – 70% of plasma protein

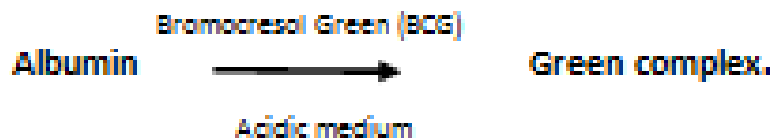- as Color's intensity increases, albumin's concentration increases

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- Albumin's reagent consists of: BCG.

• Practical Part:

- Add 3 test tubes with 1 ml reagent to each one
- Add 10 ml of the sample and of the standard to the test tubes
- Measure the absorbance , At 37 C with $\lambda_{max} = 630$
- Find the concentration by beer's law.
- standard concentration for albumin = 5 mg/dl

$\lambda_{max} = 630 \text{ nm}$

→ Let's summarize it:

We have two ways to determine the concentration of the components in the plasma:

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- Here we will try to summarize things you should remember:

Triglyceride's reagent: LPL, GK, GPO, POD

Triglyceride's λ_{max} : 500 nm

Triglyceride's standard conc.: 200 mg/dl

Glucose's reagent: GOD, POD

Glucose's λ_{max} : 500 nm

Glucose's standard conc.: 100 mg/dl

Cholesterol's reagent: Cholesterol esterase, POD.

Cholesterol's λ_{max} : 500 nm

Cholesterol's standard conc.: 200 mg/dl

Albumin's reagent: BCG (indicator)

Albumin's λ_{max} : 630 nm

Albumin's standard conc.: 5 mg/dl.

THE END

DONE BY: Hussein Alsaleh, Hisham Darwish, Dina Nasser Eddine, Lynne Abdullah, Abed Al Hadi Raed, Omar Al-Ali, Danah Alghanim and Ayham Haddad.

Questions:

1-All of the followings are used as a reagent in triglycerides experiment except:

A – LPL

B – GOD

C – GK

D – POD

2- We add to all of the following experiments a different enzymes as a catalyst to get the concentration except?

A – Triglyceride.

B – Cholesterol.

C – Glucose.

D – Albumin.

3 – The colour of Quinoneimine is directly proportional to the concentration of H_2O_2 ?

A – True.

B – False.

Answers: -

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