


Hybridoma selection



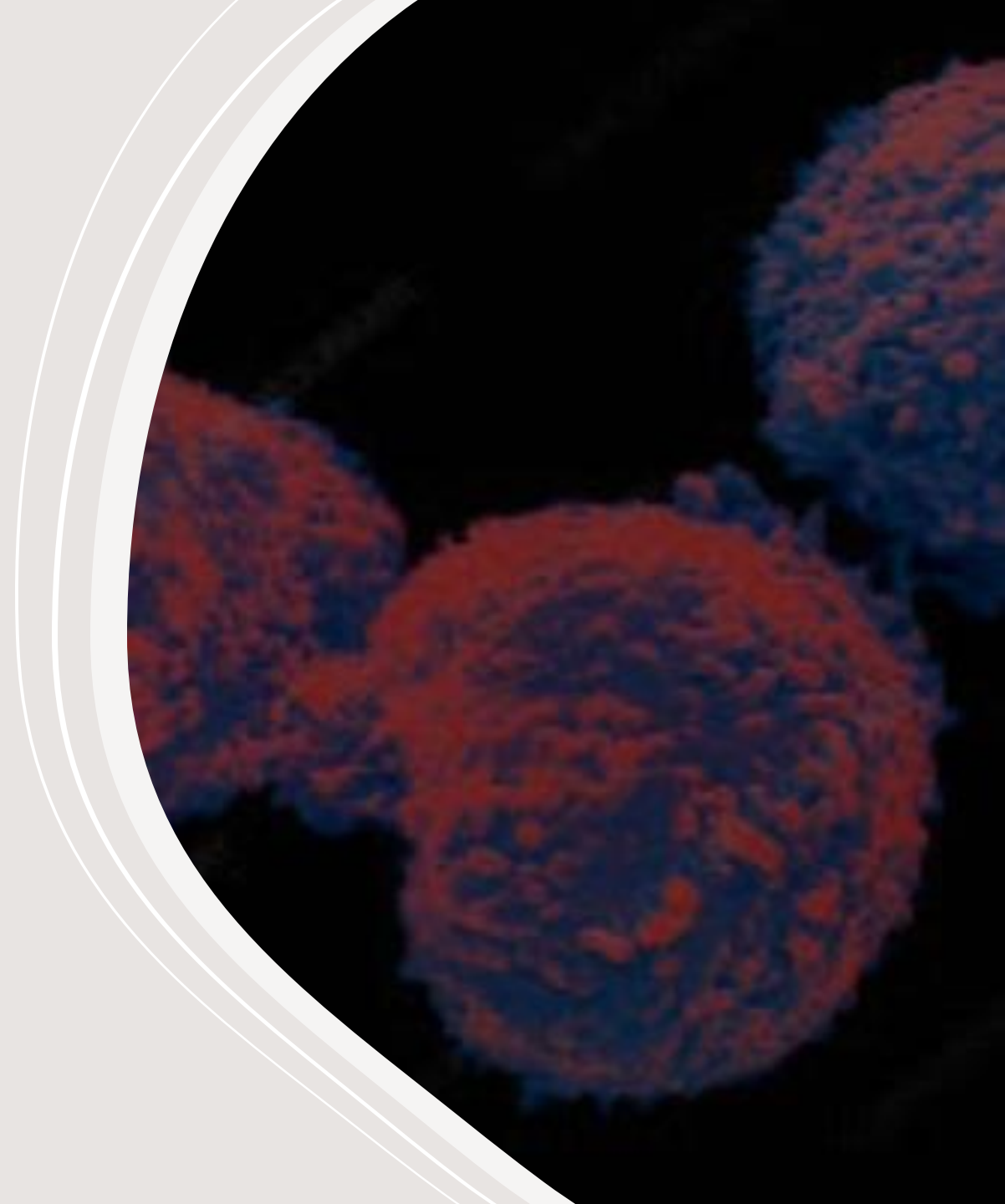

Testing of HAT medium

- Grow Myeloma cells or Hybridoma cells in HAT medium: Myeloma should die while hybridoma survive

- 
- On the day following fusion, half the old medium is removed and replaced with HAT
 - This medium change is repeated on day 2, 4, 6, 8, 11 and every 3 days afterwards
 - Within 3-4 days all myeloma cells should be dead. From the fifth day, hybrid s should start to become visible.
 - Screen for specific Ab production (11-14 day) by ELISA Select positive Ab-secreting wells.

All cells fused or unfused die (B, B-B, M, M-M) except Hybridoma cells

- The selected wells by ELISA contain the following:

A. Specific hybridoma cells: that secrete specific Abs against the injected Ag. The types of specific hybridoma cells depend on the number of epitopes in the injected Ag. If the Ag contains 5 epitopes, 5 types of hybridoma cells would be expected.

B1/M B2/M B3/M B4/M B5/M

B. Non-specific hybridoma cells: against many non-specific Ags are also present. These will be removed by the cloning method.



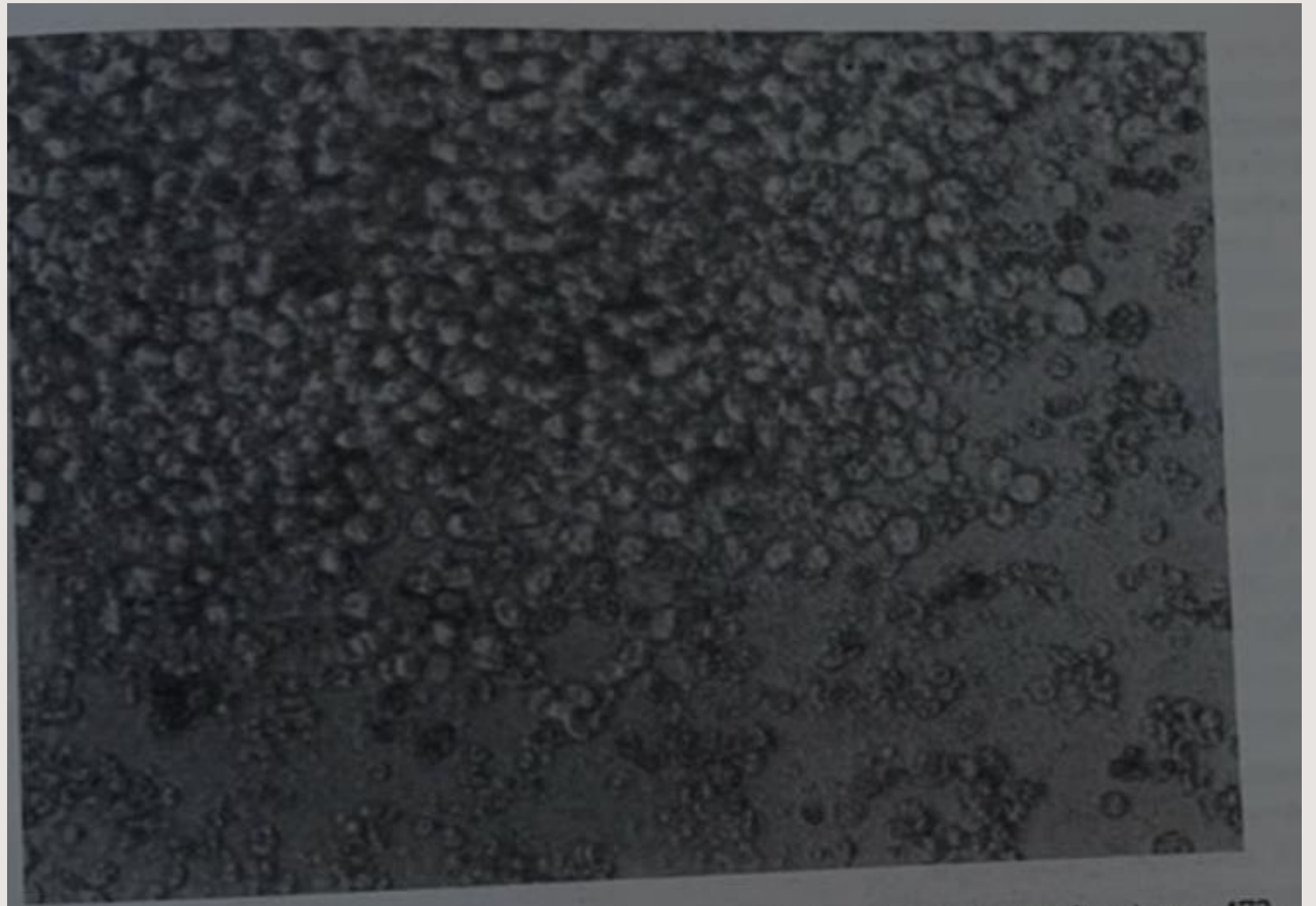


Fig. 6.1c Cell mixture after seven days in HAT medium. Magnification $\times 472$.

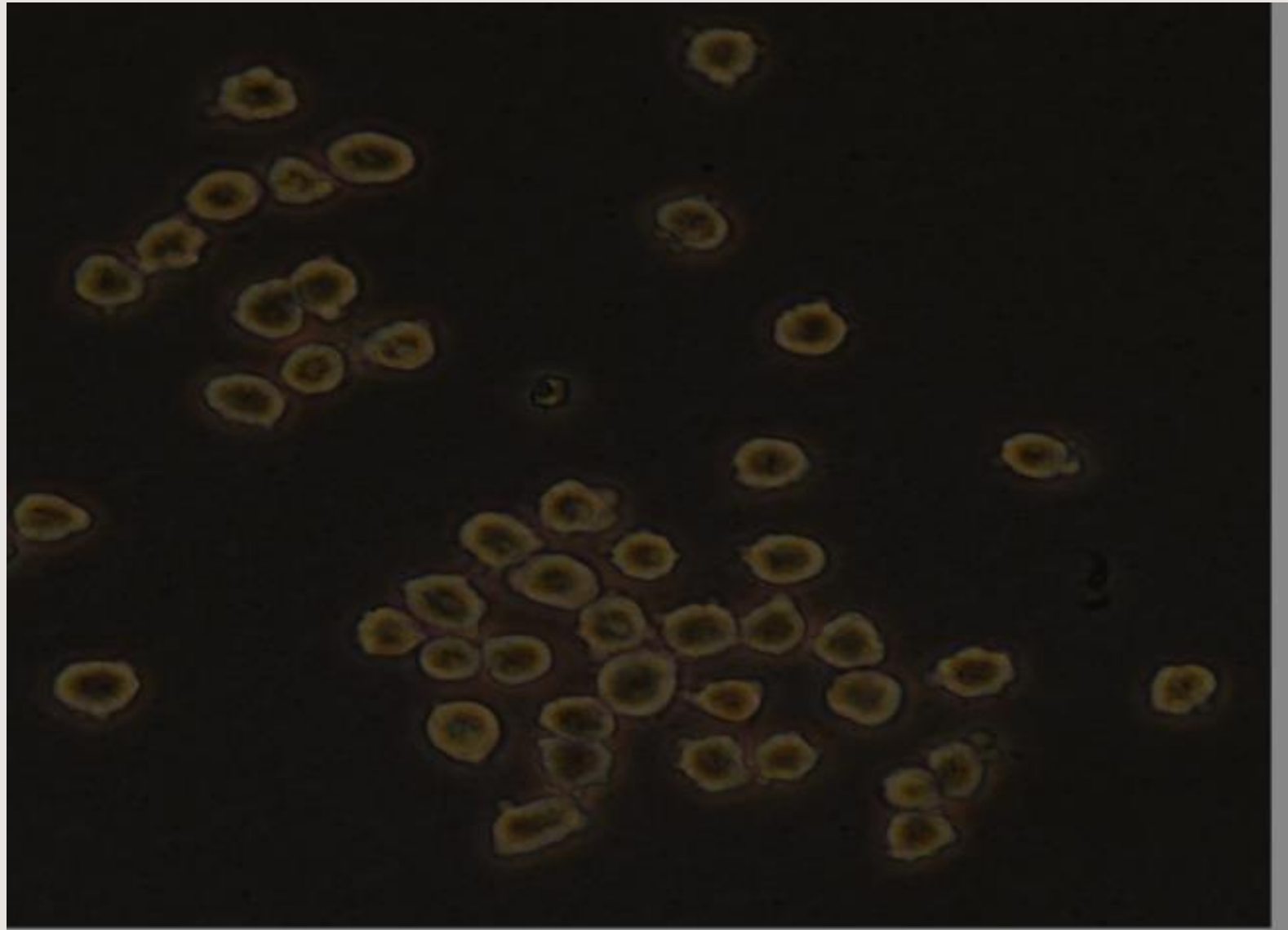

Hybridoma cells grown in tissue culture

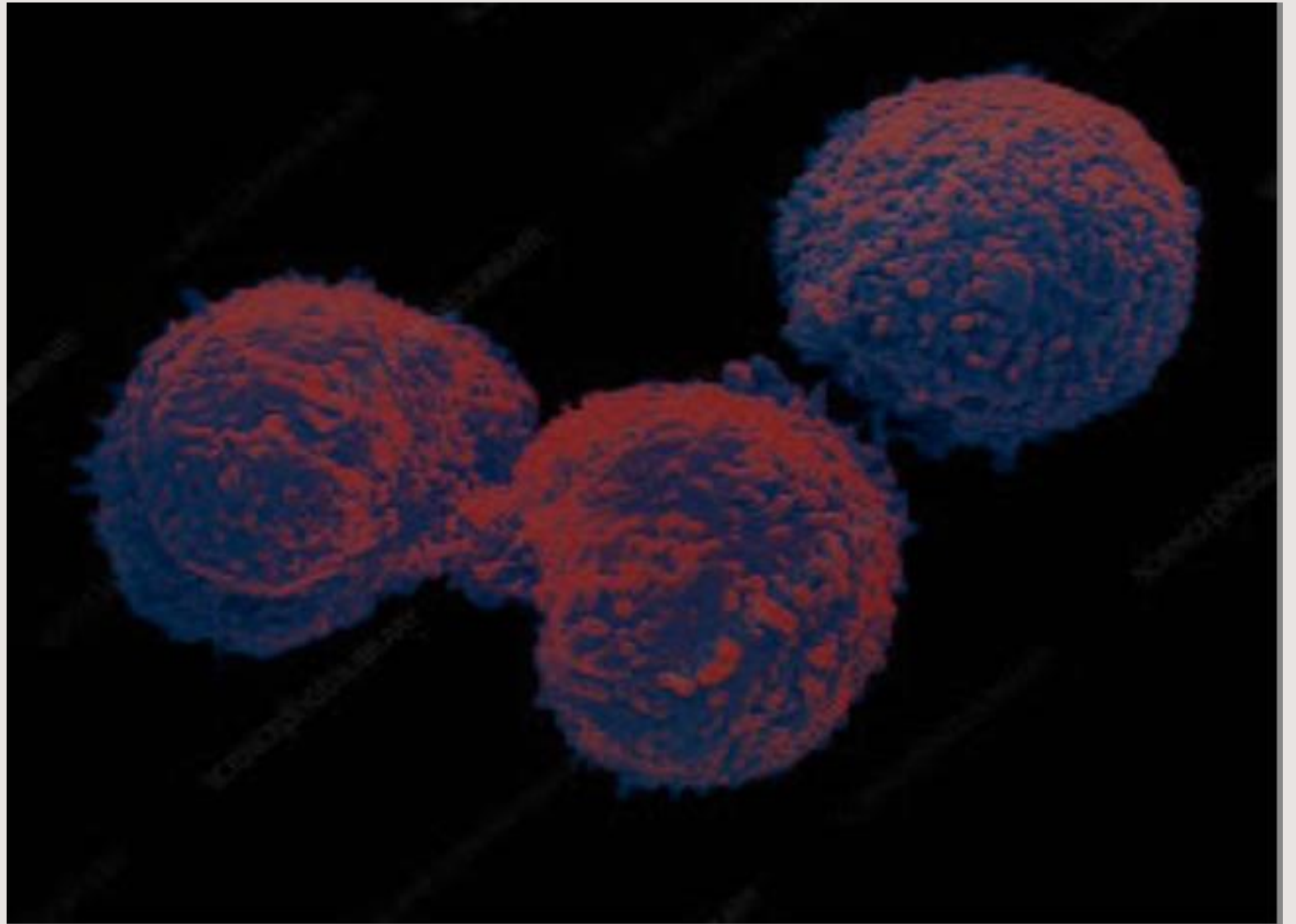

SEM of hybridoma cells